

MOBILETT XP (all products)

SP

Installation and Start-up

System (XP / Eco / Hybrid / Digital)

Installation Giraffe Stickers

© Siemens AG 2004

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Table of Contents

3

1

Prerequisites

4

General.

4

Tools required.

4

Time required

4

Parts required

4

2

Installation Procedure

5

Set of stickers.

5

Overview.

5

Location of stickers.

6

Sticker installation

7

3

Changes to the Previous Version

10

.

10

General

The application surfaces should be clean and dry prior to assembly.

The stickers should not be exposed to abnormal stretching or pulling. This could result in bubbles or a deterioration of the shape stability.

Avoid touching the adhesive side of the sticker. This could decrease the adhesive capacity of the sticker.

The best application temperature is 15-20°C, but you can still apply the stickers at temperatures as low as 10°C.

If air bubbles appear, you can puncture these with a needle and release the air using the plastic applicator or your finger.

To keep the stickers as fresh as possible, it is important to keep in mind that the stickers become optimally adhesive approximately 48 hours after assembly. Therefore, you should wait at least a couple of days before washing the stickers for the first time.

Tools required

- Masking tape

Time required

Approximately one hour for 1 person

Parts required

Installation kit, material no. 83 91 901, containing the following parts:

- Installation Instructions
- Plastic applicator
- Set of 11 stickers
- Giraffe tray

Set of stickers

Overview

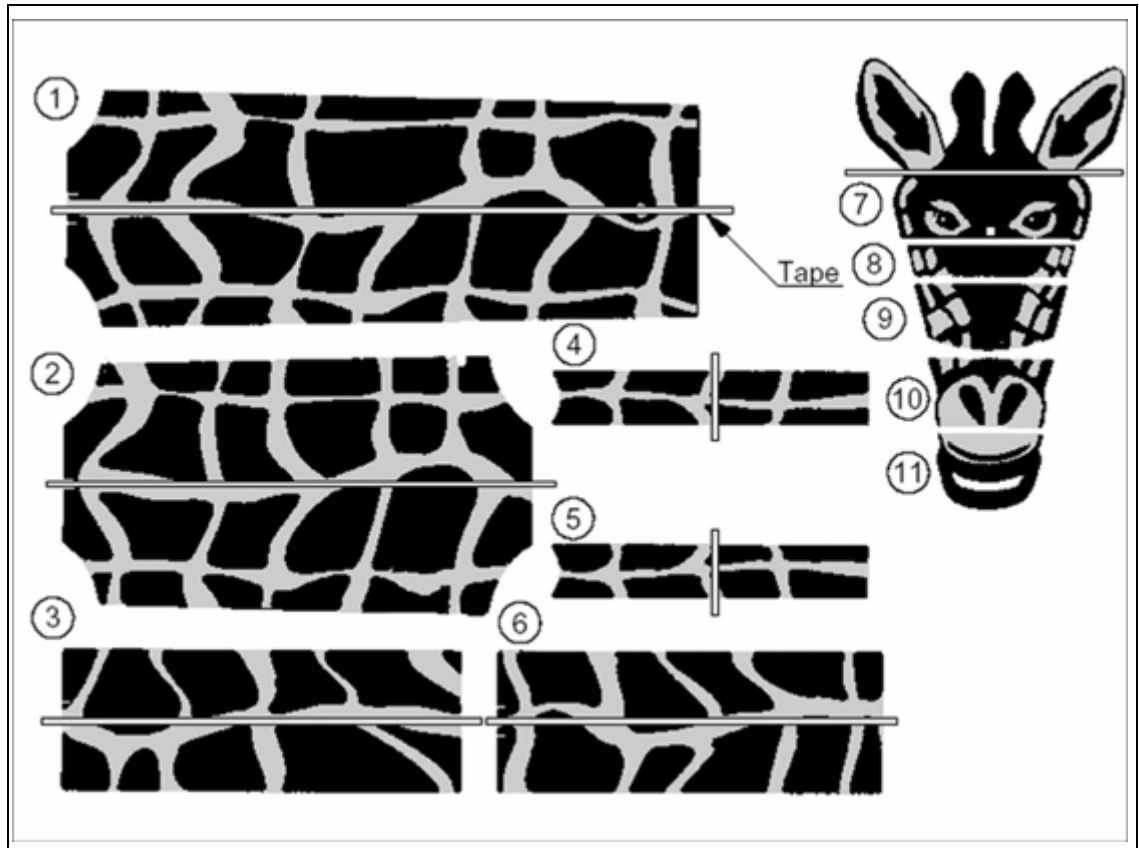


Fig. 1: Set of stickers_

The installation kit of giraffe stickers contains 11 different stickers ([Fig. 1 / p. 5](#)). The figure also shows where to place the masking tape.

Location of stickers

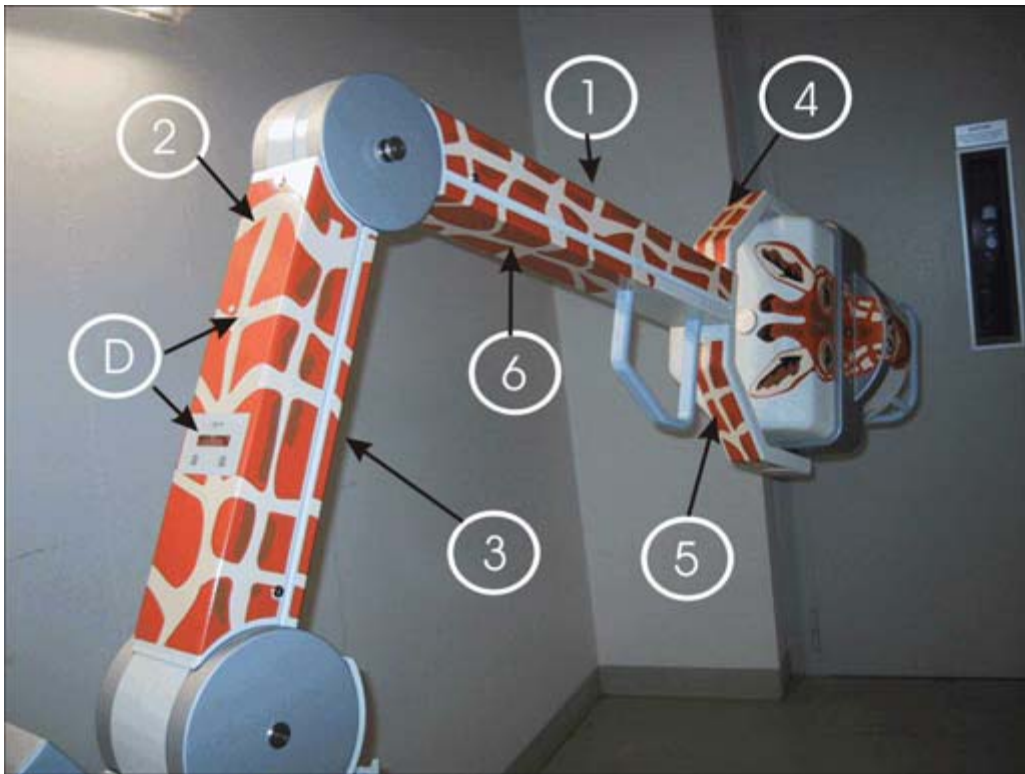


Fig. 2: Location of stickers 1_



Fig. 3: Location of stickers 2_

The installation kit of giraffe stickers contains 11 different stickers (Fig. 2 / p. 6) and (Fig. 3 / p. 6).

Sticker installation

Stickers 1-7

NOTE

When attaching sticker 2, remove the DAP display and the remote control receiver ((D/Fig. 2 / p. 6) - if these options are installed).

First use masking tape to attach the stickers to the covers as shown in (Fig. 2 / p. 6) and (Fig. 3 / p. 6). The tape will be used as a hinge. Figure (Fig. 1 / p. 5) shows where to place the tape.

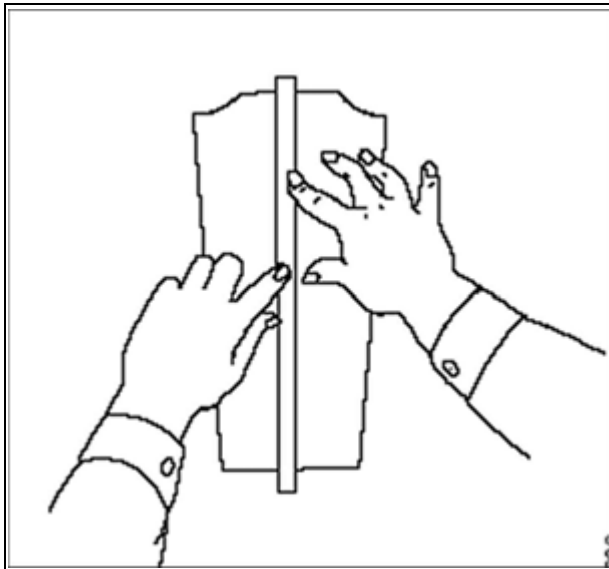


Fig. 4: Masking tape_

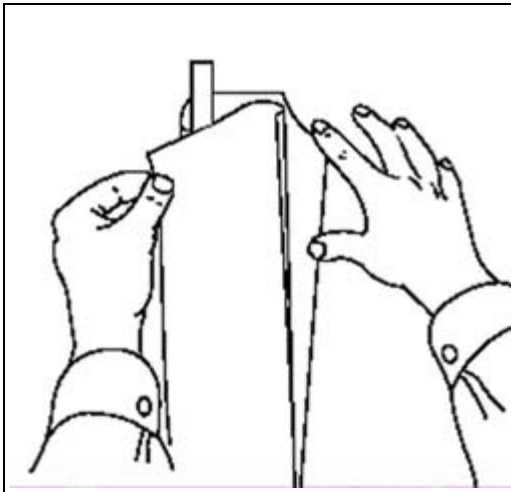


Fig. 5: Release paper_

Lift the sticker all the way to the tape hinge at an angle of 180°. Remove the release paper and carefully lower the sticker, while you paste the sticker to the cover with the other hand or with the plastic applicator ((Fig. 5 / p. 7)).

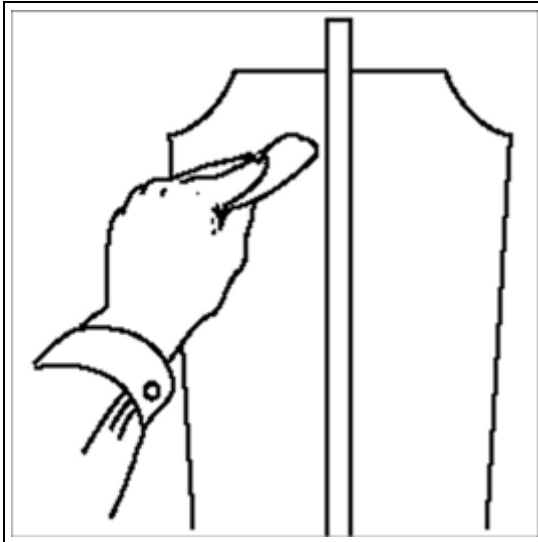


Fig. 6: Use of plastic applicator_

Use the plastic applicator from the hinge to the side. Work in overlapping motions. Remove the tape and repeat the same procedure on the other side of the sticker. If air bubbles appear, you can puncture these with a needle and release the air using the plastic applicator or your finger.

Stickers 8-11

NOTE

When attaching sticker 8, make sure that the red focus indicator is visible.

Remove the release paper and carefully paste the sticker to the covers ((Fig. 3 / p. 6)). If air bubbles occur, follow the procedure described above.

Final procedure

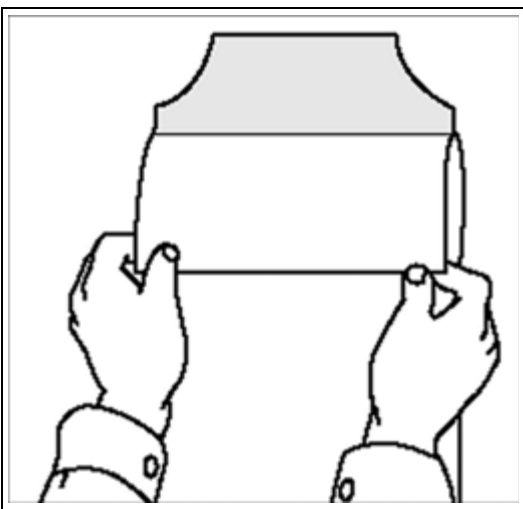


Fig. 7: Protective paper_

Remove the protective paper at an angle of 180°.

n.a., initial version